



## *Designed for Routine Weighing Applications in Your Workplace*

Offering accuracy and repeatability in essential weighing applications in laboratory, industrial and education settings, PR Analytical and Precision Balances deliver competitive performance at an economical price. Featuring RS232 connectivity for easy communication, and a backlit display and a simple interface for uncomplicated operation, the PR is perfectly designed for your workplace.

### *Standard Features Include:*

- **Basic Functionality for Routine Weighing Applications**

The Pioneer is equipped with three essential weighing modes, RS232 connectivity for data transfer and storage, and internal calibration, making it ideal for routine weighing applications.

- **Designed for Uncomplicated Operation with Easy-to-Use Display and Interface**

Equipped with an easy-to-read, bright backlit display and a simple user interface, the PR is incredibly easy to operate, with almost no training required.

- **Smart Design and Durable Construction**

The PR's small footprint saves desktop space while providing a large weighing surface. The PR is durably constructed, and features a stainless steel pan to withstand day-to-day use in the workplace.

# PR SERIES *Analytical Balances*

Standard RS232 connectivity for easy data transfer.



A large, bright backlit allows for easy viewing of results.



Durably constructed, and features a stainless steel pan to withstand day-to-day use.



Two adjustable feet that allow user to level the balance.

Four mechanical keys enable easy operation of basic application modes.

InCal Model™	PR124	PR224				
ExCal Model	PR64/E	PR124/E	PR224/E	PR223/E	PR423/E	PR523/E
Capacity (g)	62	120	220		420	520
Readability (g)	0.0001			0.001		
Repeatability (STDEV) (g)	0.0001			0.001		
Linearity (g)	0.0002			0.002		
Stabilization Time (s)	3			2		
Sensitivity Temperature Drift (PPM/K)	±3			±8	±3	
Typical Minimum Weight USP (USP K=2, U=0.10%)	200 mg			2 g		
Optimized Minimum Weight (USP, U=0.10%, K=2) SRP≤0.41d*	82 mg			0.82 g		
Units	Milligram, Gram, Kilogram, Ounce, Pound, Carat, Pennyweight, Ounce Troy, Newton, Grain					
Applications	Basic Weighing, Parts Counting, Percent Weighing					
Platform Size (diameter)	3.5 inch / 9 cm			4.7 inch / 12 cm		
Tare Range	Full range					
Power Supply	Power Input: 100 – 240V ~ 200mA 50 – 60Hz 12 – 18VA Power Output: 12 VDC 0.5A					
Assembled Dimensions (W × D × H)	208 × 320 × 309 mm					
Communication	RS232					
Operating Temperature Range	Operating conditions for ordinary lab application: +10 to 30°C (operability guaranteed between +5 and 40°C)					
Storage Temperature Range	Humidity: maximum relative humidity 80% for temperatures up to 30 °C, decreasing linearly to 50% relative humidity at 40 °C					
Net Weight	10 lb / 4.5 kg					
Shipping Weight	15.4 lb / 7 kg					
Shipping Dimensions (W × D × H)	20 × 15 × 21 inch / 507 × 387 × 531 mm					

\*The value for SRP is the standard deviation for n replicate weighing (n≥10)

# PR SERIES *Precision Balances*

Standard RS232 connectivity for easy data transfer.



A large, bright backlit allows for easy viewing of results.



Two adjustable feet that allow user to level the balance.


Four mechanical keys enable easy operation of basic application modes.

ExCal Model	PR1602/E	PR2202/E	PR4202/E	PR4201/E	PR6201/E
Capacity (g)	1600	2200	4200		6200
Readability (g)	0.01		0.1		
Repeatability (STDEV) (g)	0.01		0.1		
Linearity (g)	0.02		0.2		
Stabilization Time (s)	1				
Sensitivity Temperature Drift (PPM/K)	±6		±3		±10
Typical Minimum Weight USP (USP K=2, U=0.10%)	20 g		200 g		
Optimized Minimum Weight (g) (USP, U=0.10%, K=2) SRP≤0.41d*	8.2 g		82 g		
Units	Milligram, Gram, Kilogram, Ounce, Pound, Carat, Pennyweight, Ounce Troy, Newton, Grain				
Applications	Basic Weighing, Parts Counting, Percent Weighing				
Platform Size (diameter)	7.1 inch / 18 cm				
Tare Range	Full range				
Power Supply	Power Input: 100 – 240V ~ 200mA 50 – 60Hz 12 – 18VA Power Output: 12 VDC 0.5A				
Assembled Dimensions (W × D × H)	208 × 320 × 98 mm				
Communication	RS232				
Operating Temperature Range	Operating conditions for ordinary lab application: +10 to 30°C (operability guaranteed between +5 and 40°C)				
Storage Temperature Range	Humidity: maximum relative humidity 80% for temperatures up to 30 °C, decreasing linearly to 50% relative humidity at 40 °C				
Net Weight	7.7 lb / 3.5 kg				
Shipping Weight	11 lb / 5 kg				
Shipping Dimensions (W x D x H)	220 × 15 × 12 inch / 550 × 385 × 291 mm				

\*The value for SRP is the standard deviation for n replicate weighing (n≥10)

# PR SERIES *Analytical and Precision Balances*

## Dimensions

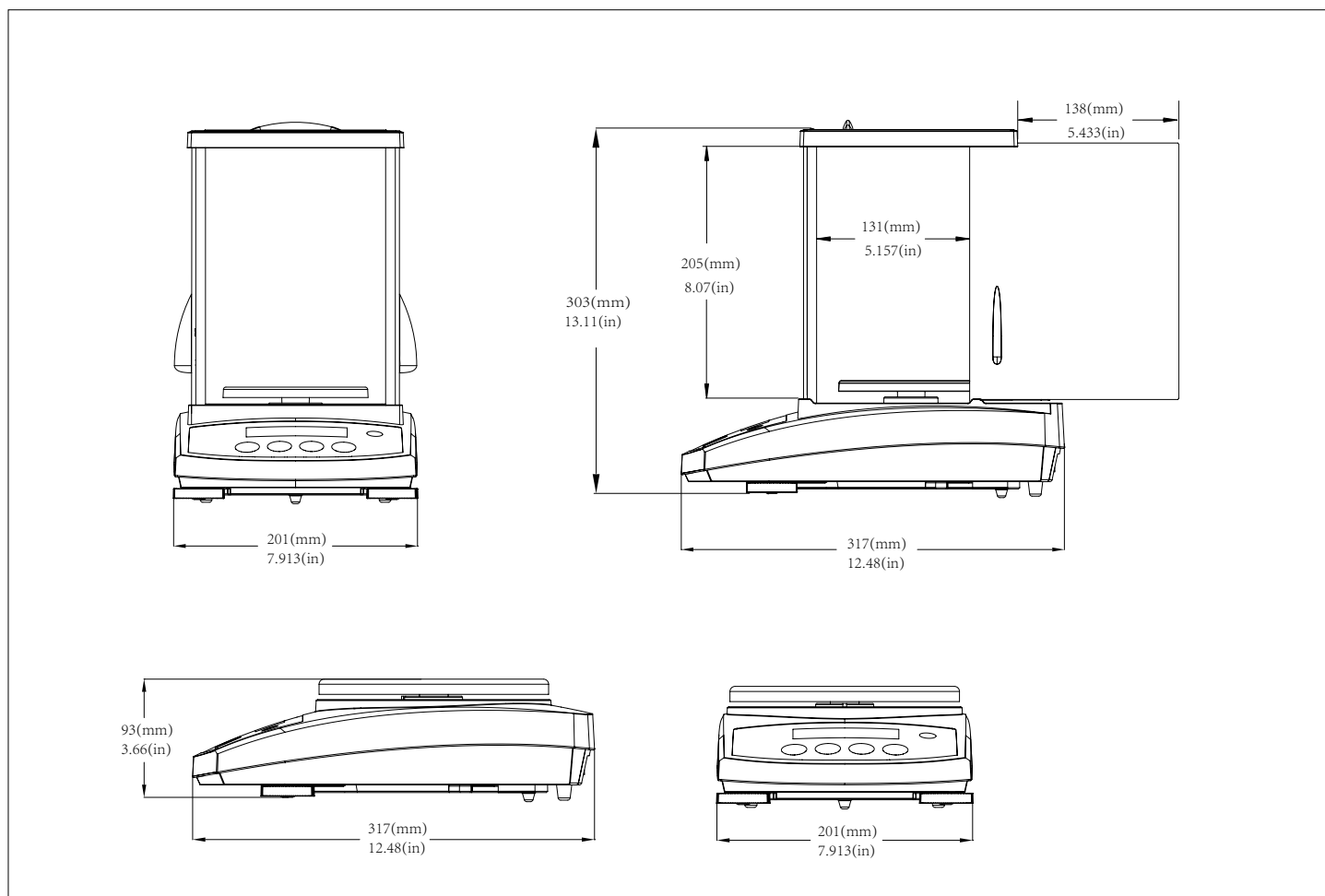

## Other Standard Features and Equipment

HIBS top housing, removable stainless steel pan, removable glass draftshield with sliding top door, integrated weigh-below-hook, security bracket, calibration lock, User-selectable environmental filters and brightness settings, auto-tare, auto-dim, user-selectable span calibration points, software lockout and reset menu, user-selectable communication settings and data print options, user-definable project and user IDs, software overload/underload indicator, stability indicator

## Compliance

- **Product Safety:** IEC/EN 61010-1; CAN/CSA C22.2 61010-1; UL 61010-1
- **Electromagnetic Compatibility:** IEC/EN 61326-1 Class B, Basic Environments; FCC Part 15 Class A; Canada ICES-003 Class A
- **Compliance Marks:** CE; CSA; RCM

## Accessories

Auxiliary Display PAD .....	80251396
Density Kit .....	80253384
Sinker for Liquid Density Determination .....	83034024
Security Device.....	80850043
RS232 Cable (25-pin) .....	80500524
Dust Cover .....	30093334
In-use Cover, (PR).....	30372547
Printer SF40A .....	30045641
Power Adapter for Balance .....	46001724

**OHAUS Asia Pacific Headquarters**  
6F, Block 7, 471 Guiping Road  
Shanghai 200233  
China

e-mail: [ChinaSales@ohaus.com](mailto:ChinaSales@ohaus.com)  
Tel: +86 21 64855408

**www.ohaus.com**

*With offices throughout Europe, Asia, and Latin America*

The management system governing the manufacture of this product is ISO 9001:2015 certified.

SAP Literature #30385150

8077264 20180209 © Copyright OHAUS Corporation

